


## Buffers, Staining and Fixing Solutions (Western Blotting)

### Ready to use Ponceau-S Stain

(For Western Blotting)

**Description:** Ponceau-S is a dye used for rapid and reversible staining of protein bands transferred on to nitrocellulose membranes, PVDF and cellulose acetate membranes. It is used to confirm transfer of protein bands from SDS-PAGE on to the membrane after electroblotting and before proceeding with immunological detection. It can also be used to mark and cut lanes on the nitrocellulose membrane before probing with different antibodies. The staining is reversible by rinsing the membrane in excess water.

**Storage:** RT

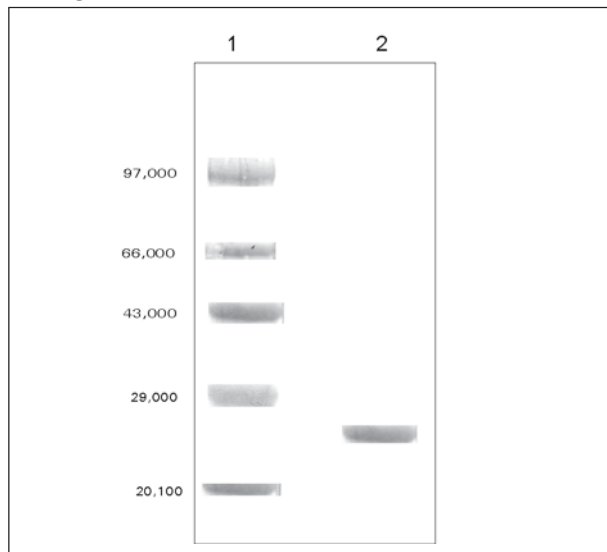

PMW-M and Protein stained with Ponceau-S and fixed with Ponceau-S fixative.

#### Ordering Information:

Product	Size	Cat #
Ready to use Ponceau-S Stain	250 ml	105373

### Ready to use Ponceau-S Fixative

(For Western Blotting)

**Description:** Protein bands transferred on to nitrocellulose membranes, PVDF and cellulose acetate membranes, stained with Ponceau-S can be fixed and preserved for reference. It is especially useful in preserving protein marker lanes on membranes when pre-stained markers are not used.

**Storage:** RT

#### Ordering Information:

Product	Size	Cat #
Ready to use Ponceau-S Fixative	250 ml	105374

### GeNei™ Western Eraser(2X)

**Description:** GeNei™ Western Eraser provides a robust and gentle method for stripping primary and secondary antibodies from blots to enable few re-probings on the same membrane.

Western Eraser buffer is ideal for use with Chemiluminescent substrates.

**Store at:** RT

#### Highlights:

- Ready and easy to use
- Compatible with commonly used Western Blotting reagents and membranes
- Saves valuable time and samples

#### Ordering Information:

Product	Size	Cat #
GeNei™ Western Eraser (2X)	500 ml	116635